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INTENDED USE

The pedal exerciser is a medical device recommended for users with mobility difficulties and is designed to restore coordination and increase muscle mass in the arms and legs.

SAFETY INSTRUCTIONS

USER: BEFORE using the Rollator, read this manual and save for future reference.

CAUTION: DO NOT attempt to adjust or operate the equipment without carefully reading all the instructions in this manual.

- **DO NOT USE OR INSTALL** the device in question for any purpose other than that expressly indicated in this manual.


- Consult your doctor or therapist to ensure proper use of the equipment.
- Keep the cardboard packaging away from sources of heat.
- The product's useful life depends on the wear and tear of its parts.
- Use with caution in the presence of children.
- All supports must be firmly supported on the surface during use to ensure the correct balance of the pedal board and the safety of the user.
- Do not use a pedal board that is damaged. There is a risk of accident!
- Unauthorised changes to the structure of the pedal board or the use of non-original parts may create dangerous situations and invalidate the warranty.
- It is recommended that the integrity of the product be checked periodically (every 3 or 4 months).

IMPORTANT SAFETY NOTICE


Do not have anyone push you while you are seated on the Rollator. This is a walking aid only and is not to be used as a transportation device.

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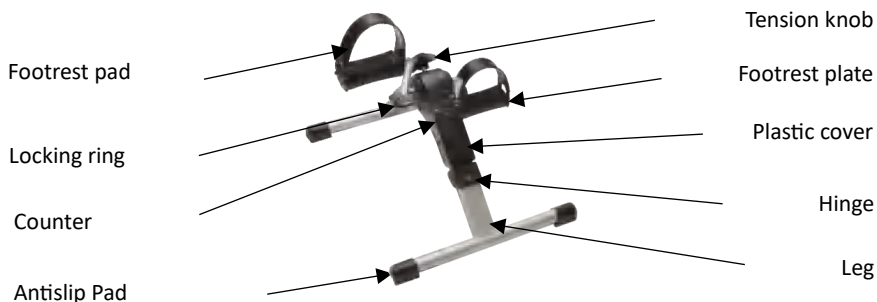
2. Pull the legs in opposite directions until the locking button clicks into the locking receiver (Figure A).

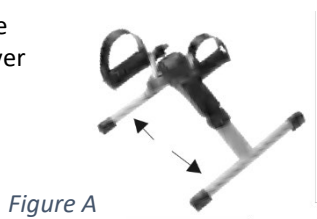3. To close, pull the locking ring and fold the leg down. Repeat on the other leg (Figure B).

Figure B

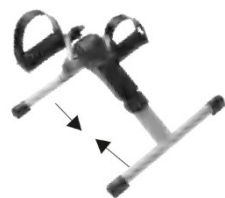4. Insert the tension knob into the top of the pedal exerciser and turn clockwise to tighten. The tension knob will increase/decrease the amount of tension on the pedal exerciser. Tighten/loosen to a level that is comfortable.

OPERATING INSTRUCTIONS

MODE (Red Button): This control allows you to select a specific function.

To start, press and hold the red button for 4 seconds. The LED will display the different functions. (If there is no activity after switching on, the unit will switch off automatically). To select a function, wait until the 4-second setting is complete. Press the red button and the arrow will move through the functions.

TIME: Determines the time the pedal exerciser has been in operation.

COUNT (CNT): Determines the number of steps during the operation of the pedal exerciser.

CALORIE (CAL): Determines the number of calories burned during the use of the pedal exerciser.

REPS PER MINUTE (R.P.M.): Automatically calculates the repetitions per minute.

SCAN: Automatically reviews the results of all functions.

Use with lower limbs

1. Place the pedal board on a flat, non-slip surface free of any objects.
2. Sit in front of the device at a distance that allows you to easily place your feet on the pedals.
3. Adjust the brake resistance according to your doctor or therapist's instructions.

NOTE: To increase the intensity of the effort, turn the brake knob clockwise; to decrease the intensity, turn it counterclockwise.

4. Operate the device by making a rotating movement.

Use with upper limbs

1. Place the pedal exerciser on a table with a flat, non-slip surface.
2. Sit in front of the device at a distance that allows you to easily operate the pedals with your hands.

3. Adjust the brake resistance according to your doctor's or therapist's instructions.

NOTE: To increase the intensity of the effort, turn the brake knob clockwise; to decrease the intensity, turn it counterclockwise.

4. Grip the pedals firmly with your hands and operate the device by making a rotating movement.

BATTERY REPLACEMENT

If the display does not light up or shows an incorrect reading, the battery needs to be replaced. The battery should only be replaced by a competent person.

MAINTAINENCE

Keep away from direct heat sources.

Use a soft, damp cloth and mild soap to clean. Dry thoroughly before use. NEVER immerse in water or other liquids. Do not clean with pressure washers or steam cleaners.

Note: Do not use acids, alkalis and/or solvents such as acetone for cleaning.

STORAGE & TRANSPORTATION CONDITION

Should be stored in clean, ventilated and dry place, rainproof, moisture-proof, and handle with care during transportation

DISPOSAL

Please dispose the product after use to comply with local regulation.

WARRANTY

This product is guaranteed for three years from the date of purchase by the user. The manufacturer will supply replacement parts for any part that has broken due to manufacturing faults. Please contact the distributor from whom you purchased the product in case you need to make a claim. Faults resulting from accidental damage, misuse, alteration or wear and tear are not covered by this guarantee. Please contact the dealer from whom you purchased the product to activate the warranty.

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